


ASSAY PERFORMANCE SUMMARY

5-acetylamino-6-formylamino-3-methyluracil (AFMU)

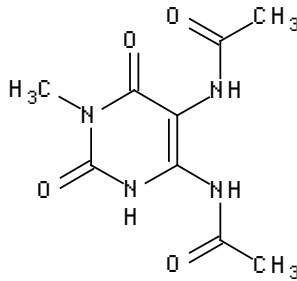

Method	LC/MS/MS
Matrix	human urine
Anti-coagulant	n/a
LLQ	100 ng/mL
ULQ	20,000 ng/mL
Sample Volume	50.0 μ L
Internal Standard	analogue
Accuracy: %Bias	
LLQ	-10.1
LQC	-6.3
MQC	-4.7
HQC	-4.7
DilQC (DF=10)	-7.2
Precision: %CV	
LLQ	14.4
LQC	8.5
MQC	5.5
HQC	7.3
DilQC (DF=10)	4.7
Stability: Minimum	
In Ice	6 hours
-20°C	244 days
No. Freeze/thaws	3