


ASSAY PERFORMANCE SUMMARY

Famotidine

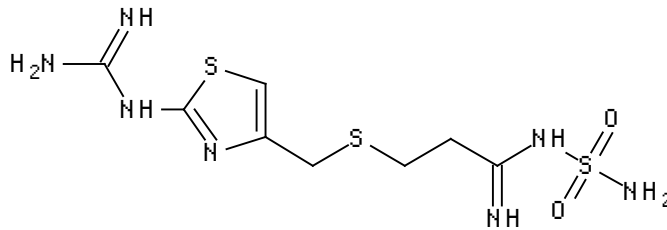

Method	LC/MS/MS
Matrix	human plasma
Anti-coagulant	K ₃ EDTA
LLQ	1.00 ng/mL
ULQ	500 ng/mL
Sample Volume	50.0 µL
Internal Standard	Stable label
Accuracy: %Bias	
LLQ	2.0
LQC	-4.3
MQC	0.7
HQC	0.0
DilQC (DF=10)	0.0
Precision: %CV	
LLQ	6.5
LQC	10.9
MQC	8.9
HQC	5.0
DilQC (DF=10)	7.6
Stability: Minimum	
Ambient	11 hours
-20°C	71 days
No. Freeze/thaws	3