


# ASSAY PERFORMANCE SUMMARY

## Fluticasone Propionate

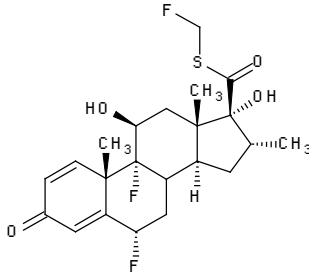

Method	LC/MS/MS
Matrix	human plasma
Anti-coagulant	K <sub>3</sub> EDTA
LLQ	5.00 pg/mL
ULQ	1000 pg/mL
Sample Volume	500 µL
Internal Standard	Stable label
Accuracy: %Bias	
LLQ	7.6
LQC	2.7
MQC	1.3
HQC	-1.8
DilQC (DF=10)	7.1
Precision: %CV	
LLQ	8.1
LQC	5.7
MQC	4.9
HQC	5.9
DilQC (DF=10)	3.3
Stability: Minimum	
Ambient	7.5 hours
-20°C	43 days
No. Freeze/thaws	4