


ASSAY PERFORMANCE SUMMARY

Methotrexate

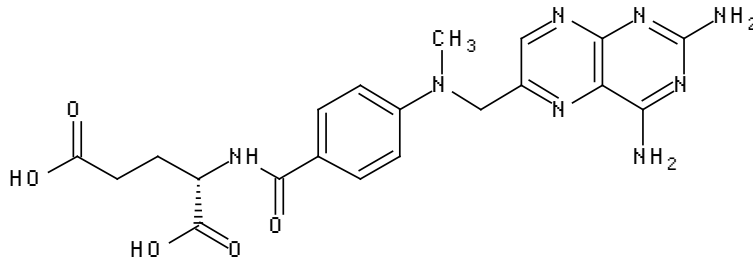

Method	LC/MS/MS
Matrix	human plasma
Anti-coagulant	NaEDTA
LLQ	1.00 ng/mL
ULQ	500 ng/mL
Sample Volume	200 µL
Internal Standard	Stable label
Accuracy: %Bias	
LLQ	4.0
LQC	2.7
MQC	0.0
HQC	-1.8
DilQC (DF=10)	0.0
Precision: %CV	
LLQ	10.5
LQC	4.5
MQC	4.3
HQC	5.0
DilQC (DF=10)	1.4
Stability: Minimum	
Ambient	6.5 hours
-20°C	527 days
No. Freeze/thaws	5