


ASSAY PERFORMANCE SUMMARY

Nelfinavir

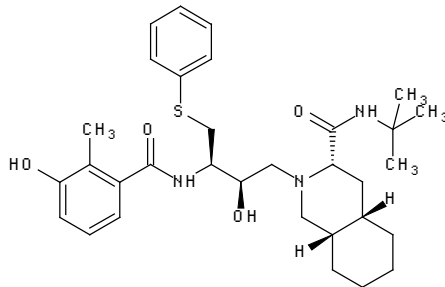

Method	LC/MS/MS
Matrix	human plasma
Anti-coagulant	K ₃ EDTA
LLQ	1.00 ng/mL
ULQ	1000 ng/mL
Sample Volume	100 µL
Internal Standard	Stable label
Accuracy: %Bias	
LLQ	6.0
LQC	9.3
MQC	5.0
HQC	3.3
DilQC (DF=10)	-0.3
Precision: %CV	
LLQ	7.3
LQC	3.7
MQC	1.4
HQC	2.5
DilQC (DF=10)	0.9
Stability: Minimum	
Ambient	6 hours
-20°C	100 days
No. Freeze/thaws	4