


ASSAY PERFORMANCE SUMMARY

Sertraline

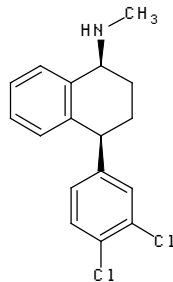

Method	LC/MS/MS
Matrix	human plasma
Anti-coagulant	Sodium heparin
LLQ	0.500 ng/mL
ULQ	250 ng/mL
Sample Volume	50.0 μ L
Internal Standard	analogue
Accuracy: %Bias	
LLQ	-2.4
LQC	-2.0
MQC	0.8
HQC	-4.0
DilQC (DF=10)	0.0
Precision: %CV	
LLQ	3.8
LQC	5.1
MQC	3.1
HQC	4.3
DilQC (DF=10)	1.4
Stability: Minimum	
Ambient	6 hours
-20°C	491 days
No. Freeze/thaws	3