


ASSAY PERFORMANCE SUMMARY

Tadalafil

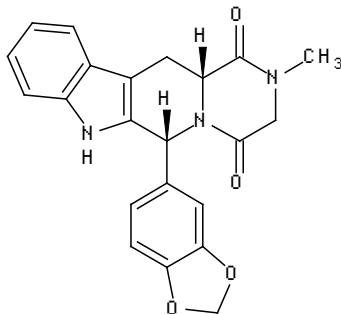

Method	LC/MS/MS
Matrix	human plasma
Anti-coagulant	K ₃ EDTA
LLQ	1.00 ng/mL
ULQ	500 ng/mL
Sample Volume	150 µL
Internal Standard	stable label
Accuracy: %Bias	
LLQ	-1.1
LQC	7.0
MQC	5.0
HQC	4.3
DilQC (DF=10)	3.0
Precision: %CV	
LLQ	11.0
LQC	3.8
MQC	2.1
HQC	2.3
DilQC (DF=10)	2.8
Stability: Minimum	
Ambient	6 hours
-20°C	137 days
No. Freeze/thaws	5